

A NEW ANGLE IN NARROW DIAMETER IMPLANTS

THE SATURNO™ NARROW
DIAMETER IMPLANT SYSTEM

Featuring an innovative, angled O-Ball
implant option which incorporates
patented, pivoting O-Ball technology.



WHAT IS YOUR ANGLE ON NARROW DIAMETER IMPLANTS?

LET ZEST ANCHORS CHANGE YOUR POINT OF VIEW!



ZEST Anchors, the trusted manufacturer of end-to-end solutions in the resilient implant-retained overdenture market is innovating again with the SATURNO™ Narrow Diameter Implant System.

SATURNO is a one-piece narrow diameter implant with an O-Ball connection featuring the patented, pivoting cap technology you are accustomed to with LOCATOR.® [This technology provides maximum resiliency while maintaining a minimized profile](#) - technology not found with the rigid connections of O-Ball mini implants on the market today.

In addition to patented pivoting technology on the restorative side, ZEST also developed a SATURNO Angled implant option to address the challenges associated with implant placement in the maxillary ridge.

Don't let patients walk out of your office - Increase your practice revenue and offer alternate treatment options like SATURNO when cost and invasive procedures inhibit treatment acceptance.

SATURNO™ — TREAT MORE PATIENTS WITH ADVANCED TECHNOLOGY

THE CASE FOR CHANGING YOUR ANGLE:

Flexibility

- Straight and Angled O-Ball options available in two cuff heights (2 & 4mm), which provide clinical flexibility
- Small O-Ball design (1.8mm) and corresponding Denture Cap require less acrylic removal thereby increasing overall denture strength
- Implant diameters available in 2.0, 2.4, or 2.9mm and in lengths of 10, 12, or 14mm enable placement in narrow ridges
- Sophisticated instrumentation for flapless surgery and drilling protocols for dense and soft bone allow you to treat more patients

Advanced O-Ball Implant Technology

- Angled option is ideal for placement in thin, off-angle ridge, providing a new treatment option
- Unique pivoting O-Ring Denture Cap enables correction of divergent implant placement and increases resiliency of the overdenture
- Proven implant design for primary stability

Practice & Patient Solutions


- Cost-conscious alternative for patients with limited finances
- Minimally invasive treatment alternative for patients who decline bone grafting




SATURNO™ VS. PREVIOUS DESIGNS

THE PIVOTAL DIFFERENCE


TRADITIONAL O-BALL MINI IMPLANT

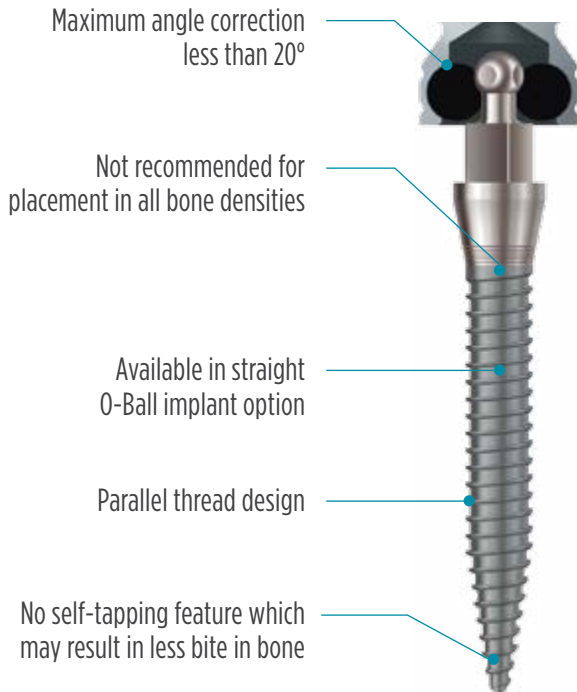

SATURNO NARROW DIAMETER IMPLANT

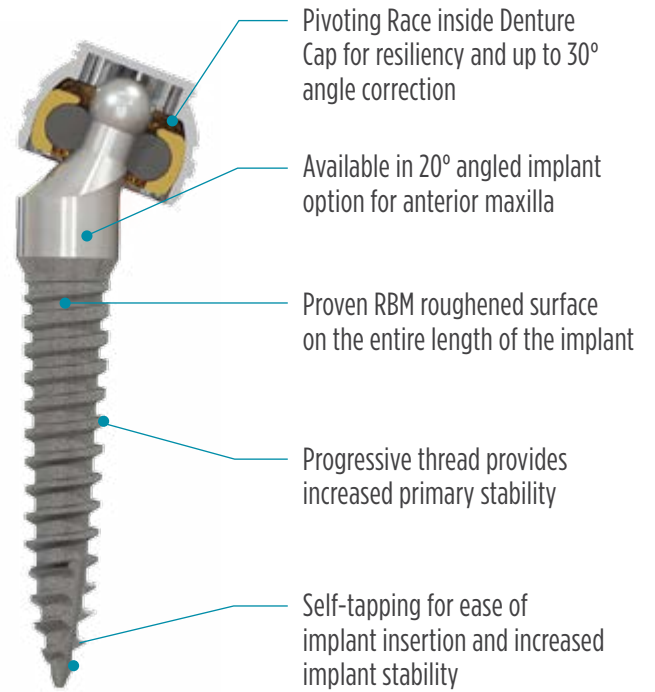

THE SATURNO O-RING ATTACHMENT OUTPERFORMS TRADITIONAL O-RINGS

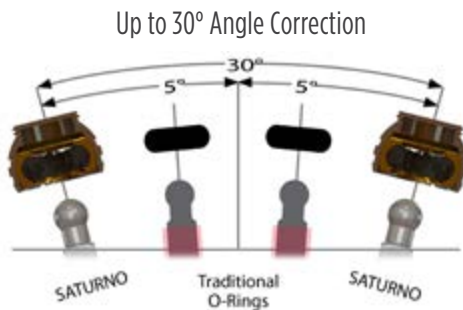

Extensive laboratory testing has demonstrated SATURNO O-Ring Attachments outperform traditional O-Rings when tested at 15° of divergence*.

*Internal data on file.

After 50,000 Test Cycles



SATURNO O-RING

ZEST Saturno O-Ring remains in excellent condition after 50,000 test cycles.




TRADITIONAL O-RING

Traditional O-Ring after 50,000 test cycles is non-functional.

THOUGHTFULLY DESIGNED INSTRUMENTATION

SATURNO™ NARROW DIAMETER IMPLANT SYSTEM

Feel confident using SATURNO Instrumentation. SATURNO osteotomy preparation instruments have easily identifiable laser etched depth markings for clear visibility when drilling and placing implants. A drill stop set is included for an additional level of depth control. SATURNO instruments were designed to enable predictable osteotomy preparation, especially with a flapless protocol. In addition, every SATURNO Implant is packaged with a convenient driver mount for seamless delivery to the implant osteotomy.

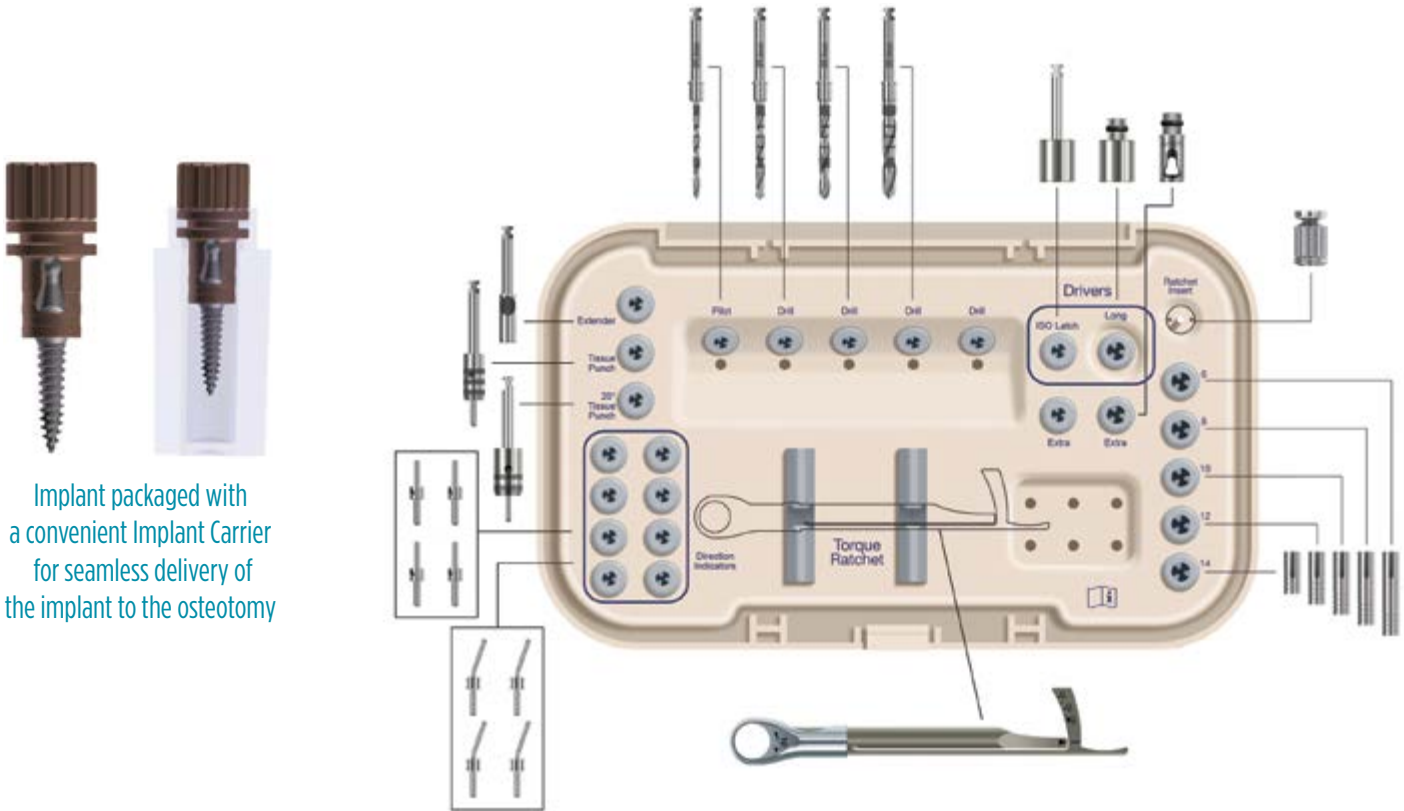

Implant packaged with a convenient Implant Carrier for seamless delivery of the implant to the osteotomy

SATURNO ALL-IN-ONE PACKAGING



ANGLED OR STRAIGHT

Your choice of 2.0, 2.4 or 2.9mm Diameters

Your choice of 2 or 4mm Cuff Heights



Restorative Component Vial



Micro O-Ring attachment assembly (O-Ring cap, O-Ring race & O-Ring)



Block-out Spacer



SATURNO™
NARROW DIAMETER IMPLANT SYSTEM



CHANGE YOUR
PERSPECTIVE ON
NARROW DIAMETER
IMPLANTS BY OFFERING
your edentulous patients the SATURNO™
Narrow Diameter Implant System today!



ZEST ANCHORS
2061 WINERIDGE PLACE
ESCONDIDO, CA 92029
USA

TEL: 1.800.262.2310
FAX: 760.743.7975
EMAIL: SALES@ZESTANCHORS.COM
WWW.ZESTANCHORS.COM

CE
0086

DISTRIBUTED BY: