




**LOCATOR<sup>®</sup>**  
OVERDENTURE IMPLANT SYSTEM

---

THE ONLY IMPLANT ENGINEERED FOR THE WORLD'S  
MOST TRUSTED OVERDENTURE ATTACHMENT

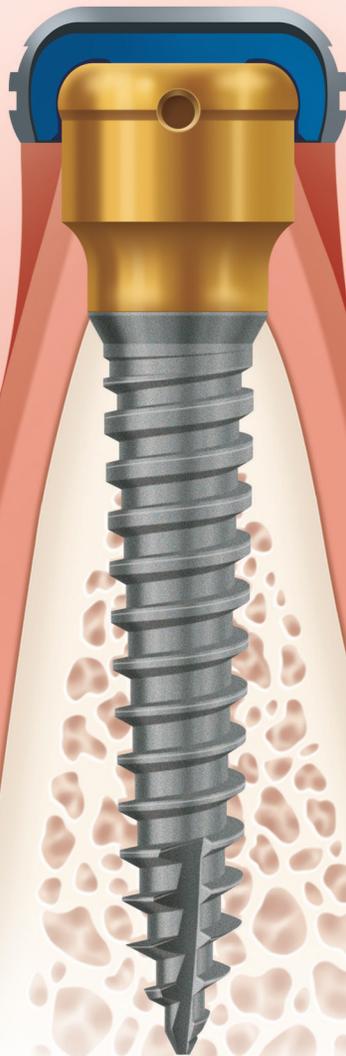

INDICATED FOR FULLY OR PARTIALLY EDENTULOUS PATIENTS  
— NARROW AND STANDARD RIDGE —


---

Zest has answered clinician requests and patient demands for a cost-conscious solution for securing overdentures. The LOCATOR® Overdenture Implant System brings quality, ease and affordability together so you can bring confidence back to more patients' lives.

## NARROW RIDGE LOCATOR IMPLANTS

Our narrow ridge implants make overdenture treatment possible for patients with severe resorption and very narrow ridges. The typical treatment plan for these edentulous patients would be bone grafting, requiring a long process with additional surgeries and added expense.

For these reasons, patients often decline dental implant treatment, reducing their overall quality of life. The narrow ridge LOCATOR implant can be placed through the tissue without flapping, allowing you to offer patients a cost-effective solution with implant retained overdentures.

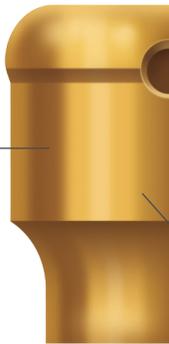

1

The Zest LOCATOR Abutment's dramatically reduced vertical height provides patient comfort when the overdenture is removed, as well as increased overdenture strength as compared to O-Ball attachments

2

Aggressive thread design that widens at the coronal portion providing primary stability when immediate loading may be indicated

3

Proven RBM roughened surface on entire length of the implant

4

Self-tapping design for ease of implant insertion and increased implant stability



LOCATOR® Narrow and Standard Ridge Implants, built for the best-in-class LOCATOR Abutment, are the clear choice to offer a less invasive, predictable and durable implant-retained overdenture option for patients.



5

The Zest LOCATOR Male self-aligns and pivots inside the Denture Cap providing a genuine resilient connection that holds-up to patient mastication forces while providing attachment durability

6

All-Inclusive packaging with industry leading Zest LOCATOR Abutments for 2.5 and 4mm cuff heights. For implant placement flexibility, individual replacement LOCATOR Abutments are also available in 3.0, 5.0 and 6mm cuff heights

7

Available in 2.4, 2.9 and 3.4mm diameters and 10, 12 or 14mm lengths for placement in all ridges

## STANDARD RIDGE LOCATOR IMPLANTS

NEW to our offering, and now available, these two-part standard ridge implants allow you to reach even more patients with the LOCATOR system. Standard ridge implants are typically used in younger patients that have more bone volume available.

### ALL-INCLUSIVE PACKAGE CONTAINS

- 1 Implant
- 1 LOCATOR Abutment
- 1 Processing Pack


Each Processing Pack has what you need to select retention levels and address draw correction; improving ease of denture placement and removal



### STANDARD MALES




### EXTENDED RANGE MALES

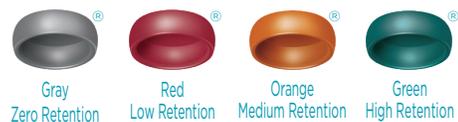



**LOCATOR®**  
OVERDENTURE IMPLANT SYSTEM

**THE LOCATOR OVERDENTURE IMPLANT SYSTEM SURGICAL KIT**  
CONTAINS PREPARATION INSTRUMENTS WITH DEPTH MARKINGS AND A DRILL STOP SET FOR SIMPLICITY AND CONTROL DURING OSTEOTOMY PREPARATION.

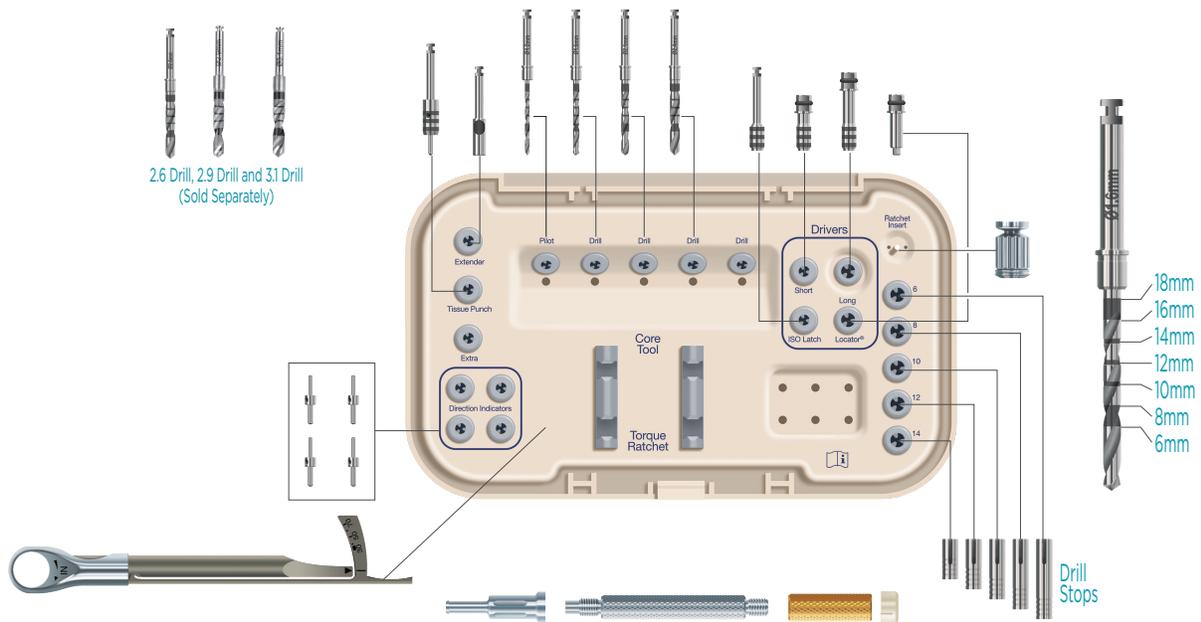

**ZD ZEST DENTAL SOLUTIONS®**

**ZEST DENTAL SOLUTIONS**  
2875 LOKER AVENUE EAST  
CARLSBAD, CA 92010  
USA

TEL: 1.855.866.LODI (5934)  
FAX: 760.743.7975  
EMAIL: SALES@ZESTDENT.COM  
[WWW.ZESTDENT.COM](http://WWW.ZESTDENT.COM)

DISTRIBUTED BY: